

NOTES ON UPGRADING TO PROCAL / PROCAL-TRACK V4x FROM V2x



**IF YOU ARE RUNNING VERSION 2x OF PROCAL OR PROCAL-TRACK
PLEASE USE THE V3x DISC TO UPGRADE FIRST, THEN RUN THE V4x UPGRADE**



V2x TO V3x FIRST, THEN V3x TO V4x



PROCAL / PROCAL-TRACK V2x TO V4x UPGRADE GUIDE

UPGRADING YOUR INSTALLATION FROM PROCAL V3x

To update an existing installation of ProCal to the latest V4x version, use the **UPGRADE** installation option on the Transmille Software CD.

This upgrade should be performed on each computer the programs are installed on.




**IF YOU ARE RUNNING VERSION 2x OF PROCAL OR PROCAL-TRACK
PLEASE CONTACT TRANSMILLE BEFORE PROCEEDING**




The **UPGRADE** will preserve any data already on the computer, only adding new files required for the new functions available.

STEP 1 : Click OK to proceed

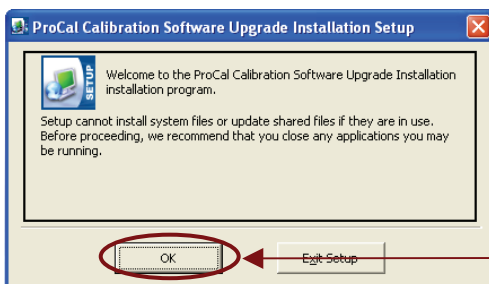

Click to Proceed


STEP 2 : Click INSTALL to proceed

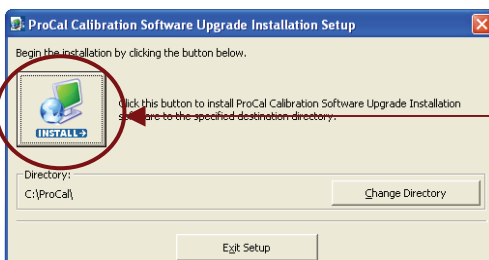

Click to Proceed




STEP 3 : Click CONTINUE to proceed

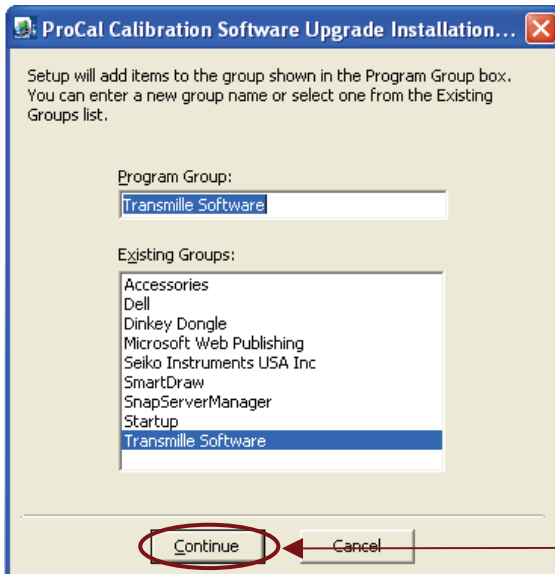

Click to Proceed

STEP 4 : Click INSTALL to proceed

ProCal will now be upgraded, displaying a progress screen as shown below :

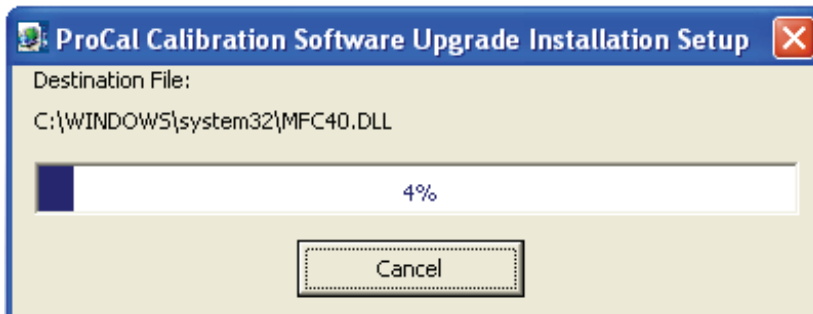

On completion of installation, the following screen will be displayed :



STEP 5 : Upgrade ProCal-Track (If Installed)



**Repeat steps 1-4
for the ProCal-track option if this program is installed.**

UPDATING YOUR FILES TO V4x



IMPORTANT STEP :: After the upgrade installer has completed, you must run the 'Updater' program – follow the steps below to complete :

STEP 1 : Start the UPDATER program


Select from the Windows Start menu :

All Programs -> Transmille Software -> ProCal-ProCal-Track Database Updater V3x to V4x




STEP 2 : Set Paths & Run Updater

Check program paths are correct (program auto detects installations, so default should be correct)
Click 'Update' button to proceed

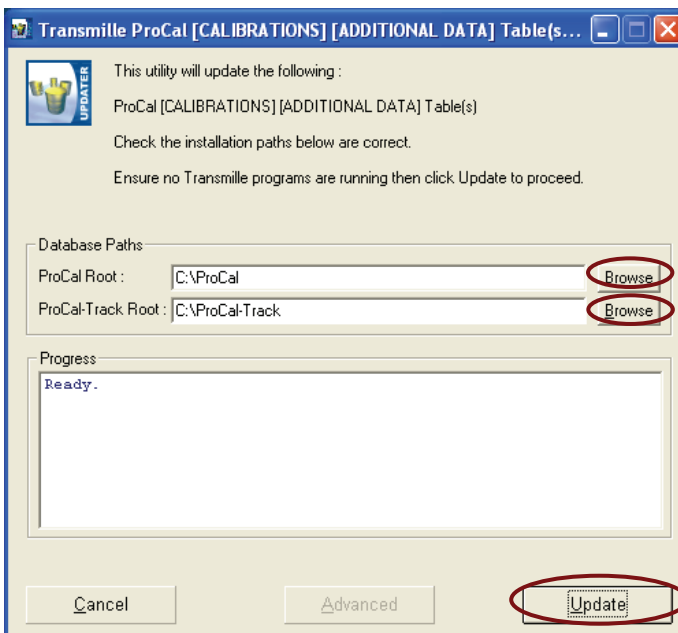


**Set program paths
or leave set as default
(DEFAULT RECOMMENDED)**

Click to start update

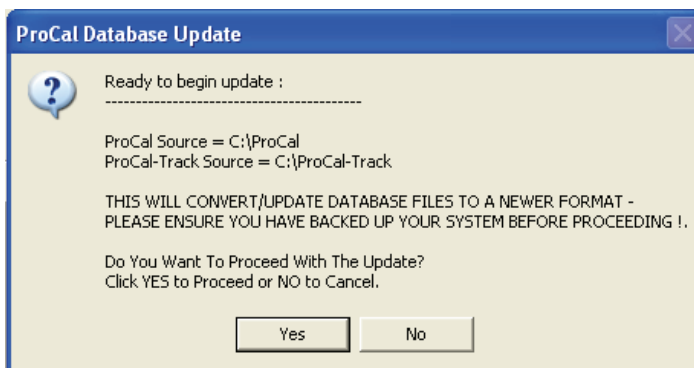

Click **YES** to proceed with update - the updater will proceed to update database files, and display progress on the screen above.

Once completed the message below is displayed :

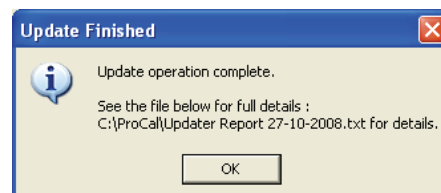

Click  to exit the updater after it has finished – update is now complete.