

NOTES ON UPGRADING TO PROCAL / PROCAL-TRACK V4x



This guide contains information on how to upgrade ProCal / ProCal-Track from V3x to V4x.



IF YOU ARE RUNNING **VERSION 2x** OF PROCAL OR PROCAL-TRACK
PLEASE CONTACT TRANSMILLE BEFORE PROCEEDING

BACKING UP PROCAL



**BEFORE PROCEEDING, PLEASE TAKE THE TIME TO PERFORM A FULL BACKUP OF
YOUR CURRENT PROCAL AND PROCAL-TRACK INSTALLATIONS.**

These installations are relatively small and the complete program folders should be backed up :

C:\ProCal

C:\ProCal-Track

PLUS ANY NETWORK / SERVER FOLDERS

These backups can be made to CD, a USB memory stick, Network Backup Storage or similar.

Once the backup is performed, carry on to the next page and follow the steps to upgrade your programs.

UPGRADING YOUR INSTALLATION FROM PROCAL V3x

To update an existing installation of ProCal to the latest V4x version, use the **UPGRADE** installation option on the Transmille Software CD.

This upgrade should be performed on each computer the programs are installed on.



**IF YOU ARE RUNNING VERSION 2x OF PROCAL OR PROCAL-TRACK
PLEASE CONTACT TRANSMILLE BEFORE PROCEEDING**




**UPGRADE OPTION
PROCAL**

**UPGRADE OPTION
PROCAL-TRACK**




The **UPGRADE** option will preserve any data already on the computer, only adding new files required for the new functions available.

STEP 1 : Click OK to proceed

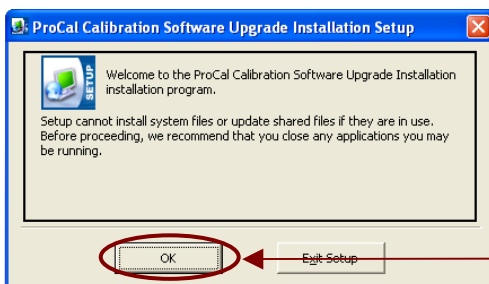


Click to Proceed

STEP 2 : Click INSTALL to proceed

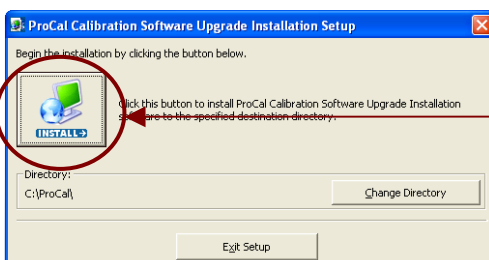

Click to Proceed


STEP 3 : Click CONTINUE to proceed

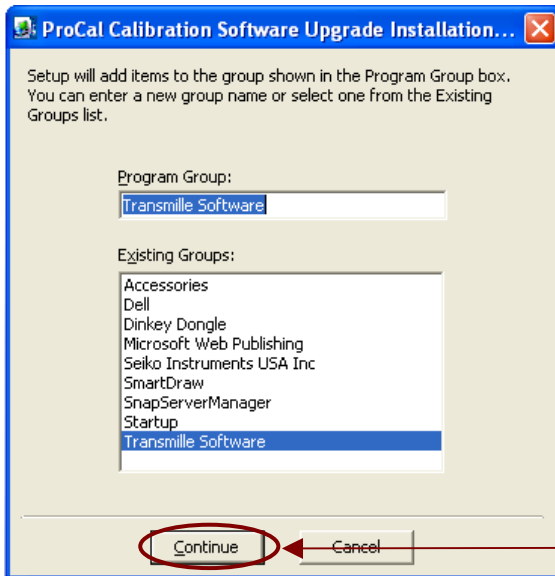


Click to Proceed

STEP 4 : Click INSTALL to proceed

ProCal will now be upgraded, displaying a progress screen as shown below :

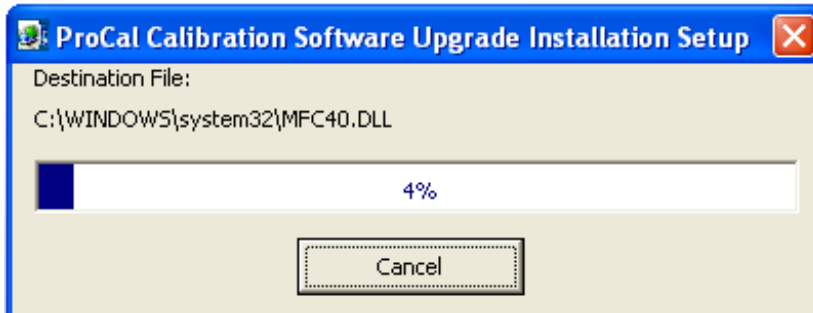

On completion of installation, the following screen will be displayed :

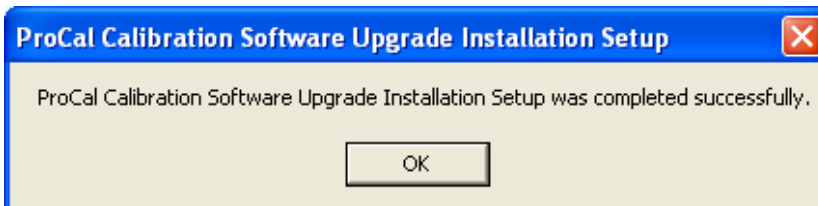

STEP 5 : Upgrade ProCal-Track (If Installed)



**Repeat steps 1-4
for the ProCal-track option if this program is installed.**

UPDATING YOUR FILES TO V4x



IMPORTANT STEP :: After the upgrade installer has completed, you must run the 'Updater' program – follow the steps below to complete :

STEP 1 : Start the UPDATER program


Select from the Windows Start menu :

All Programs -> Transmille Software -> ProCal-ProCal-Track Database Updater V3x to V4x




STEP 2 : Set Paths & Run Updater

Check program paths are correct (program auto detects installations, so default should be correct)
Click 'Update' button to proceed

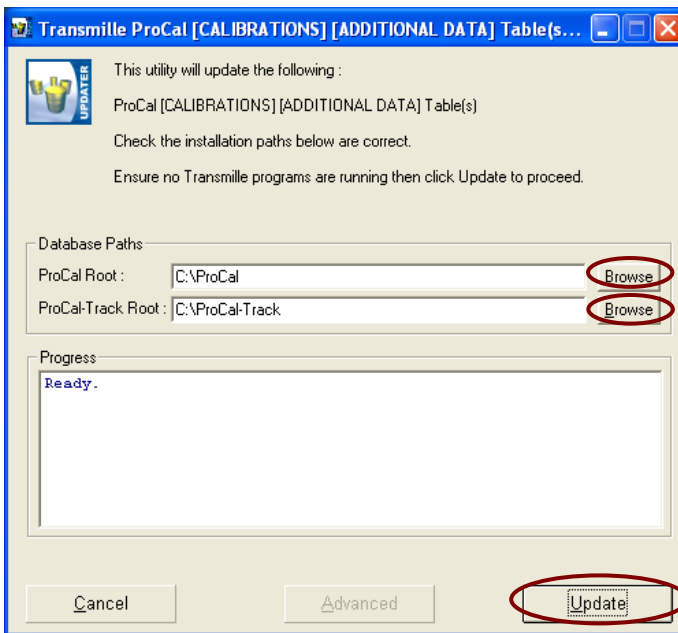


Set program paths or leave set as default (DEFAULT RECOMMENDED)

Click to start update

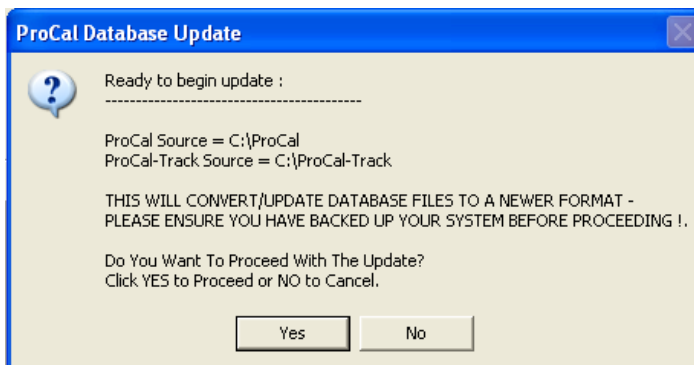

Click **YES** to proceed with update - the updater will proceed to update database files, and display progress on the screen above.

Once completed the message below is displayed :

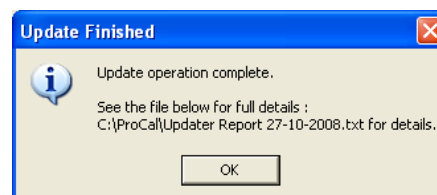

Click to exit the updater after it has finished – update is now complete.